

DF-550

Gastric-Thoracic Drainage Pump



www.drfriend.com



DOCTOR'S FRIEND

CE 2460 ISO 13485

DF-550

Gastric-Thoracic Drainage Pump



The DOCTOR'S FRIEND gastric-thoracic drainage suction pump is engineered to provide suction function by alternating the expansion and contraction of air within a cylinder at regular intervals. Simplicity of design to enhance safety and hygiene.



- The machine offers two functions - continuous and intermittent types with timer which helps the user on post-operative care.
- The digital timer adopts micro computer programming design and large displayed numbers for easy reading and simple operation, while solid-state controls help to assure mild suction function.




- The safety device is equipped within the bottle. When the bottle is full, the floating control will activate and block the adapter to prevent overflow.
- The vacuum bottle is made of impact resistant materials and graduated with 1800 c.c. capacity.

Fill in water up to water sealed line before operation.

360° moving caster with locking brake features easy movement and low center of gravity brings high safety.



Whole ABS plastic injection made outer shell and also features with carrying handle which is easy to lift and transport.

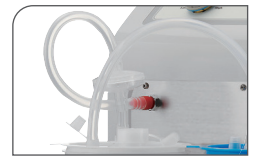

The filter may prevent the pump from pollution and helps to prolong its life.



Depositing box creates extra space for storage.



The dual-fuse device design.



Standard Accessory

Suction tube 180cm	1 pc
PC adaptor	1 pc
Filter (L)	2 pcs

Specification

Voltage	AC 240V
Vacuum range	0~1000 mm/H ₂ O
Flow rate	10 L/min
Safety device	Floating control
Capacity	1800 c.c. + 1800 c.c.
N.W. / G.W.	9/13 kg
Machine size	L36 x W33 x H60 cm
Packing size	L50 x W48 x H66 cm

90-006

Vacuum Pump



DOCTOR'S FRIEND medical instrument co., ltd.

No.28, Ln. 13, Sec. 3, Zhongshan Rd., Tanzi Dist., Taichung City 42756, Taiwan
TEL : +886-425347777 FAX : +886-425338669 Email : medical.df@doctorsfriend.com.tw

